

Libra 3405M / 3412M

The Libra 3405M / 3412M are developed with global shutter CMOS sensors, offering both high speed and high resolution. They feature a wide spectral range of 350 nm–1100 nm, providing superior performance for multi-channel fluorescence imaging. Their compact design also facilitates seamless integration into instrument systems.



Key Features

Benefits

Global Shutter	High-speed, artifact-free imaging, ideal for high-throughput slide scanning.
350 nm–1100 nm wide spectrum	High sensitivity in low light and NIR, ideal for multi-channel fluorescence.
3.4 μm Pixel Size	Optimized for < 40X optical systems, enables higher resolution for cellular and tissue imaging.
10 GigE Interface	Provides higher data throughput and more stable transmission than USB 3.0. ^[1]
Compact Design	Facilitates seamless integration into instrument systems.

Typical Applications


- Pathology Slide Scanning
- Microscopy Imaging
- Industrial Inspection

Noted Examples


[1] The maximum full-resolution speed is 163 fps for the Libra 3405M and 98 fps for the Libra 3412M, about three times faster than the USB 3.0 version.

Model	Interface	Frame Rate (fps)
Libra 3412	GigE	98 fps
	USB 3.0	28 fps
Libra 3405	GigE	164 fps
	USB 3.0	65 fps

Quantum Efficiency

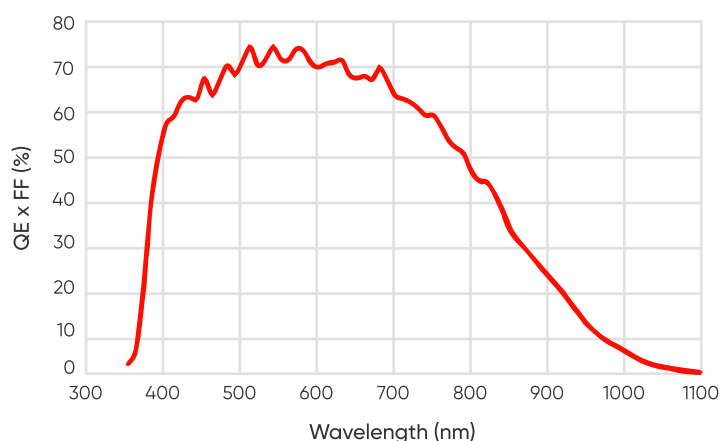

Dimensions (Unit: mm)

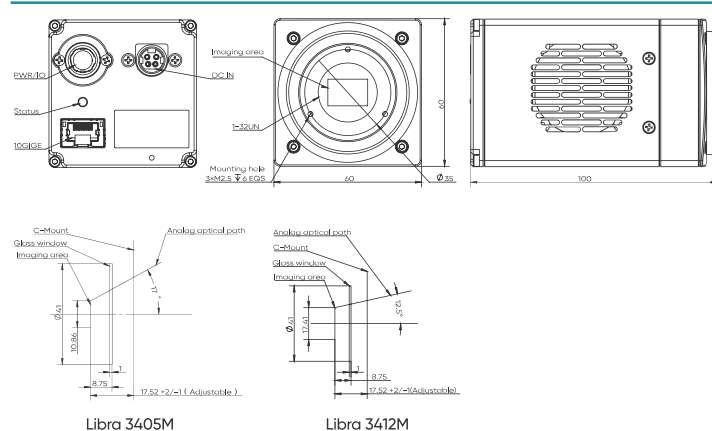

Specifications

Product Name	Libra 3405M			Libra 3412M		
Model	L3405M			L3412M		
Sensor Type	CMOS					
Sensor Model	Gpixel GMAX 3405			Gpixel GMAX 3412		
Chrome	Mono					
Array Diagonal	10.9 mm (2/3")			17.4 mm (1.1")		
Effective Area	8.3 mm x 7.0 mm			14.0mm x 10.5mm		
Pixel Size	3.4 μm x 3.4 μm					
Resolution	2448 (H) x 2048 (V)			4096 (H) x 3072 (V)		
Peak QE	75%@540 nm; 47%@800 nm					
Gain Mode	HighCapacity, Balanced, Sensitive					
Full Well Capacity	12bit: High Capacity 8.9 Ke-, Balanced 4.2 Ke-, Sensitive 0.48 Ke-			12bit: High Capacity 9 Ke-, Balanced 4.5 Ke-, Sensitive 0.7 Ke-		
Frame Rate	12 bit	10 bit	8 bit	12 bit	10 bit	8 bit
	100 fps	163 fps	163 fps	62 fps	65 fps	98 fps
Readout Noise	12 bit (Median): 3.7 e-@High Capacity, 2.3 e-@Balanced, 1.4 e-@Sensitive			12 bit (Median): 3.6 e-@High Capacity, 2.4 e-@Balanced, 1.6 e-@Sensitive		
Shutter Mode	Global					
Exposure Time	1 μs ~10 s					
Image Correction	DPC					
ROI	Support					
Binning (FPGA)	1x1 , 2x2 , 4x4					
Cooling Method	Air cooling					
Cooling Temp.	10°C@Ambient 25°C					
Dark Current	0.5 e-/p/s@10°C					
Trigger Mode	Hardware, Software					
Trigger Output	High, Low, Exposure out, Readout, Trigger Ready					
Trigger Interface	Hirose-12-Pin					
Data Interface	10 GigE					
Bit Depth	High Depth (12 bit), Standard (10 bit), Speed (8 bit)					
Optical Interface	C Mount					
Power Supply	12 V / 8 A					
Power Cons.	30 W					
Dimensions	60 mm (H) x 60 mm (W) x 100 mm (L)					
Weight	~489 g					
Software	Mosaic V3, Sample Pro, Micro-Manager 2.0					
SDK	C / C++ / C# / Python					
Operating System	Windows, Linux					
Environment	Working: Temp. 0°C~40°C, HUM 10%~85%; Storage: Temp. -10°C~60°C, HUM 0%~85%					

*Specifications in this manual are subject to changes without prior notice.



Follow us

+86 (591) 2805 5076 ext.818

www.tucsen.com

support@tucsen.com